


# PHOSPHOLIPASE A<sub>2</sub> [PLA<sub>2</sub>]

from *Streptomyces violaceoruber*  
(Phosphatidylcholine 2-acylhydrolase, EC 3.1.1.4)

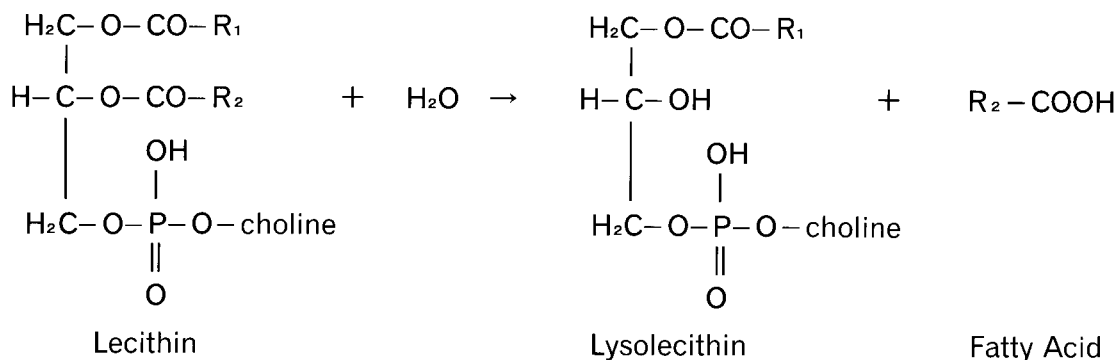

## Preparation and Specification

Appearance : White to light brownish amorphous powder, lyophilized  
Specific activity : More than 1 U/mg solid

## Properties

Substrate specificity : See Table 1  
Molecular weight : 15 kDa (Sephadex G-100)  
Isoelectric point : pH 7.51  
Michaelis constant : Lecithin  $5.0 \times 10^{-3}$  M  
Optimum pH : 7.3-8.3 Figure 1  
pH stability : 6.0-9.5 (55 °C, 10 min) Figure 2  
Thermal stability : Stable at 50 °C and below  
(pH 8.0, 10 min, + CaCl<sub>2</sub>) Figure 3  
Stabilizer : Ca<sup>2+</sup>  
Activators : Ca<sup>2+</sup>, Non-ionic detergents  
Inhibitor : EDTA

## Applications for Diagnostic Test

This enzyme is useful for enzymatic determination of **lecithin** when coupled with lysophospholipase (T-32), glycerophosphorylcholine phosphodiesterase (T-33) and choline oxidase (T-05).




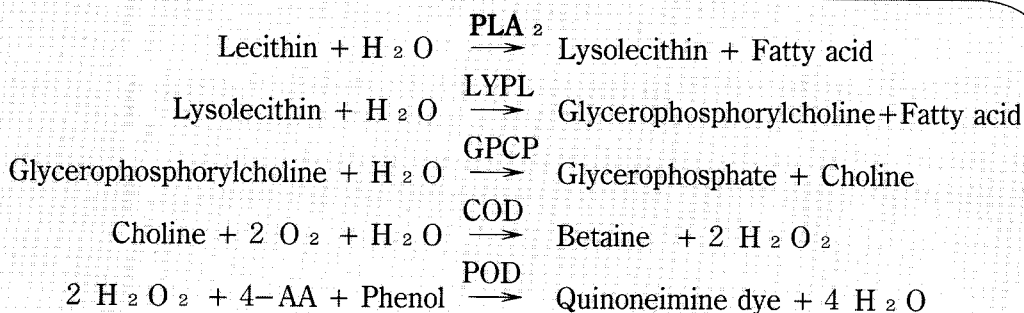


Table 1. Substrate specificity


Substrate	Relative activity (%)
Phosphatidylcholine	100
Phosphatidylethanolamine	30
Phosphatidic Acid	9.8
Lysophosphatidylcholine	0
Triolein	0
Trilaurin	0

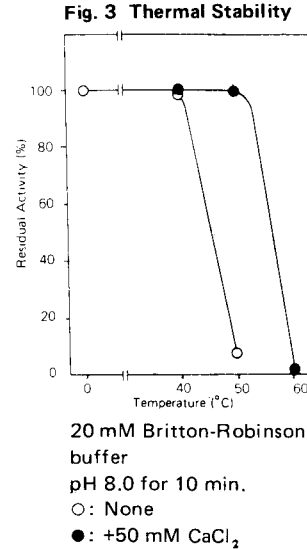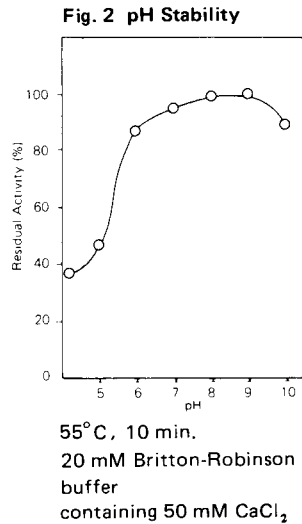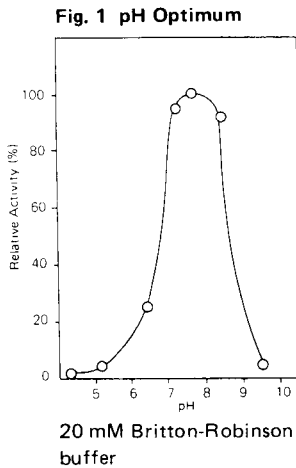

## Assay

### Principle

The assay is based on the increase in absorbance at 500 nm as the formation of quinoneimine dye proceeds in the following reactions:

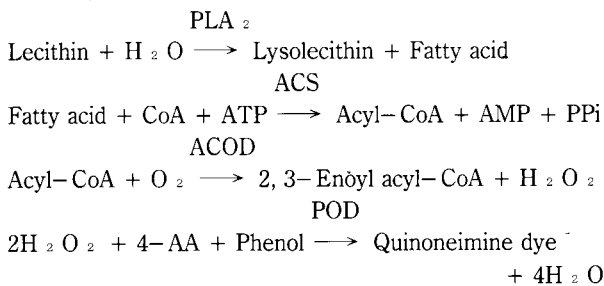

ACS: Acyl-CoA synthetase

ACOD: Acyl-CoA oxidase

ATP: Adenosine triphosphate

### Unit definition

One unit is defined as the amount of enzyme which liberates 1 μ mole of fatty acid per minute at 37°C under the conditions specified in the assay procedure.

### Reagents

1. Reaction mixture for the first reaction	
0.2 M Tris-HCl buffer pH 8.0	0.20 ml
20 mM DPPC solution <sup>1)</sup>	0.05 ml
1 M CaCl <sub>2</sub> solution	0.025 ml
0.2 M ATP solution pH 7.0	0.01 ml
0.1 M CoA solution pH 6.0 <sup>2)</sup>	0.02 ml
10.6 U/ml ACS solution <sup>3)</sup>	0.05 ml

Distilled water 0.145 ml

1): 20 mM DPPC solution

Dissolve 147 mg of DPPC with 10 ml of 2% (W/V) Triton X-100 solution.

2): 0.1 M CoA solution pH 6.0

Dissolve 7.85 g (purity calculation) of CoA with 80 ml of distilled water and adjust pH to 6.0 (25°C) with 4 N NaOH and add distilled water to make a total of 100 ml.

3): 10.6 U/ml ACS solution

Dissolve 106 U of ACS with 10 ml of Tris-HCl buffer pH 8.0.

DDPC: Dipalmytoylphosphatidylcholine

2. Reaction mixture for the second reaction

0.5 M PIPES-NaOH buffer pH 7.5 0.04 ml

0.3% (W/V) 4-AA solution 0.10 ml

0.2% Phenol solution 0.10 ml

50 U/ml POD solution <sup>4)</sup> 0.045 ml

300 U/ml ACOD solution <sup>5)</sup> 0.05 ml

10% (W/V) Triton X-100 solution 0.01 ml

0.2 M ATP solution pH 7.0 0.05 ml

10 mM FAD solution 0.005 ml

Distilled water 0.10 ml

4): 50 U/ml POD solution

Dissolve 500 U (PPU) of POD with 10 ml of distilled water.

5): 300 U/ml ACOD solution

Dissolve 3,000 U of ACOD with 10 ml of PIPES-NaOH buffer pH 7.5.

FAD: Flavine adenine dinucleotide

3. 20 mM NEM  
NEM: N-ethylmaleimide
4. Reaction stopper  
0.5 % (W/V) SDS solution containing 0.1 M EDTA  
SDS: Sodium dodecyl sulfate  
EDTA: Ethylenediamine tetraacetic acid
5. Enzyme dilution buffer  
10 mM Tris-HCl buffer pH 8.0 containing 0.05 % (W/V) BSA.
6. Reagents:  
DPPC: Nippon Fine Chemical Co., Ltd.  
Coenzyme A (CoA · Li<sub>3</sub>): Kohjin Co.  
ATP (2Na · 3H<sub>2</sub>O): Kyowa Hakko Co., Ltd.  
PIPES [Piparazine-N, N'-bis (2-ethanesulphonic acid)]: Dojindo Laboratories # 345-02225  
Triton X-100: Yoneyama Yakuhin Kogyo Co., Ltd.  
# 34021600  
NEM: Wako Pure Chemical Industries, Ltd.  
# 058-02061  
ACS: Asahi Kasei Corporation # T-16  
ACOD: Asahi Kasei Corporation # T-17  
SDS: NACALAI TESQUE, INC. # 316-06  
4-AA: NACALAI TESQUE, INC. Special grade  
# 01907  
POD: Sigma Chemical Co. Type II #P-8250  
EDTA(2Na · 2H<sub>2</sub>O): Daiichi Pure Chemicals Co., Ltd.  
Special grade

#### ■ Enzyme solution

Weigh about 20 mg of test sample exactly and add enzyme dilution buffer to make a total of 20 ml. Dilute it with enzyme dilution buffer to adjust the concentration as required.

#### ■ Procedure

- Pipette accurately 0.50 ml of reaction mixture for the first reaction into a small test tube and preincubate at 37 °C .
- After 5 min, add exactly 50 μ l of enzyme sample and mix to start the first reaction at 37 °C .  
※ In the case of a test blank, add 50 μ l of enzyme dilution buffer in place of enzyme solution at this point.

- After 10 min, add 0.5 ml of 20 mM NEM and add 0.5 ml of reaction mixture for the second reaction after 15 seconds later and mix to start the second reaction at 37 °C .
- At 5 min after starting the reaction, add 1.5 ml of the reaction stopper to stop the reaction.
- Measure the absorbance at 500 nm.

$$\begin{aligned} \text{Absorbance sample} &: A_s \\ \text{blank} &: A_b \\ \Delta A &= (A_s - A_b) \leq 0.200 \text{ Abs} \end{aligned}$$

#### ■ Calculation

$$\text{Activity (U/mg)} = \frac{\Delta A/10}{12.0 \times 1/2} \times \frac{3.05}{0.05} \times \frac{1}{X}$$

12.0 : millimolar extinction coefficient of quinoneimine dye at 500 nm (cm<sup>2</sup> / μ mole)

1/2 : a multiplier derived from the fact that 2 mole of H<sub>2</sub>O<sub>2</sub> produces 1 mole of quinoneimine dye

10 : reaction time (min)

3.05 : final volume (ml)

0.05 : volume of enzyme solution (ml)

X : concentration of the sample in enzyme solution (mg/ml)

#### Storage

Storage at -20 °C in the presence of a desiccant is recommended.

#### References

- Scandella, C. J. and Kornberg, A. C. (1971) *Biochemistry*, **10**, 4447-4456.
- Nishijima, M., Akamatsu, Y. and Nojima, S. (1974) *J. Biol. Chem.*, **249**, 5658-5667.
- Waite, M. and Sisson, P. (1971) *Biochemistry*, **10**, 2377-2383.
- Kannagi, R. and Koizumi, K. (1979) *Biochem. Biophys. Acta*, **556**, 423-433.
- Van Dam-Mieras, M. C. E., Slotboom, A. J., Pieterse, W. A. and de Haas, G. H. (1975) *Biochemistry*, **14**, 5387-5393.

## P L A<sub>2</sub> 活性測定法 (Japanese)

### I. 試薬液

#### 1. 第一反応試薬混合液

0.2M トリス-HCl 緩衝液 pH8.0	0.20 ml
20mM DPPC 溶液 <sup>1)</sup>	0.05 ml
1M 塩化カルシウム溶液	0.025 ml
0.2M ATP 溶液 pH7.0	0.01 ml
0.1M CoA 溶液 pH6.0 <sup>2)</sup>	0.02 ml
10.6U/ml ACS 溶液 <sup>3)</sup>	0.05 ml
精製水	0.145 ml

#### 1): 20mM DPPC 溶液

DPPC 147mg を 2%(W/V) トリトン X-100 溶液 10ml で溶解する。

#### 2): 0.1M CoA 溶液 pH6.0

CoA 7.85g (純度換算) を精製水 80ml に溶解した後、4N NaOH で pH6.0(25 °C) に調整し、精製水で全容 100ml とする。

#### 3): 10.6U/ml ACS 溶液

ACS 106 単位 (U) を 10mM トリス-HCl 緩衝液 pH8.0 10ml で溶解する。

#### 2. 第二反応試薬混合液

0.5M PIPES-NaOH 緩衝液 pH7.5	0.04 ml
0.3%(W/V) 4-AA 溶液	0.10 ml
0.2%(W/V) フェノール液	0.10 ml
50U/ml POD 溶液 <sup>4)</sup>	0.045 ml
300U/ml ACOD 溶液 <sup>5)</sup>	0.05 ml
10%(W/V) トリトン X-100 溶液	0.01 ml

0.2M ATP 溶液 pH7.0 0.05 ml  
 10mM FAD 溶液 0.005 ml  
 精製水 0.10 ml

4): 50U/ml POD 溶液  
 POD 500 単位 (PPU) を精製水 10ml で溶解する。

5): 300U/ml ACOD 溶液  
 ACOD 3,000 単位 (U) を 10mM PIPES-NaOH 緩衝液 pH7.5 10ml で溶解する。

3. 20mM NEM 溶液

4. 反応停止液  
 0.1M EDTA を含む 0.5%(W/V) SDS 溶液

5. 酵素溶解希釈用液  
 0.05%(W/V) BSA を含む 10mM トリス-HCl 緩衝液 pH8.0

6. 試薬  
 DPPC(ジパルミトイルフォスファチジルコリン)  
 : 日本精化製

CoA(コエンザイム・Li<sub>3</sub>): 興人製  
 ATP(アデノシン三リン酸・2Na・3H<sub>2</sub>O):  
 協和発酵製

PIPES [ピペラジン-N,N'-ビス(2-エタンスルホン酸)]: 同仁化学製 #345-02225

トリトン X-100: 米山薬品工業製 #34021600

FAD(フラビンアデニンジヌクレオチド・2Na):  
 協和発酵製

NEM(N-エチルマレイミド): 和光純薬工業製  
 特級 #058-02061

ACS(アシル-CoA シンセターゼ):  
 旭化成製 #T-16

ACOD(アシル-CoA 酸化酵素): 旭化成製  
 #T-17

SDS(ドデシル硫酸ナトリウム):  
 ナカライテスク社製 #316-06

4-AA: ナカライテスク社製 特級 #01907

POD: シグマ社製 Type II #P-8250

EDTA(エチレンジアミン四酢酸・2Na・2H<sub>2</sub>O):  
 第一化学薬品製 特級

## II. 酵素試料液

検品約 20mg を精密に量り、酵素溶解希釈用液で全容 20ml とする。

その液を酵素溶解希釈用液で適宜希釈する。

## III. 測定操作法

1. 小試験管に第一反応試薬混合液 0.50ml を正確に注入し、37℃で予備加温する。
2. 5分経過後、酵素試料液 50 μl を正確に加えて混和し、37℃で第一反応を開始する。  
 ※盲検は酵素試料液の代わりに酵素溶解希釈緩衝液 50 μl を加える。
3. 10分経過後、20mM NEM 溶液 0.50ml を加えて混和し、15秒後に第二反応試薬混合液 0.50ml を加えて混和し、37℃で第二反応を開始する。
4. 5分経過後、反応停止液 1.50ml を加えて混和し、反応を停止する。
5. 500nm における吸光度を測定する。  
 求められた吸光度を試料液は As、盲検液は Ab とする。

$$\Delta A = (A_s - A_b) \leq 0.200 \text{ Abs.}$$

## IV. 計算

$$\text{活性 (U/mg)} = \frac{\Delta A / 10}{12.0 \times 1/2} \times \frac{3.05}{0.05} \times \frac{1}{X}$$

12.0: キノンイミン色素の 500nm におけるミリモル分子吸光係数 (cm<sup>2</sup>/μmole)

1/2: H<sub>2</sub>O<sub>2</sub> 2モルからキノンイミン色素 1モルが生成することによる係数

10: 反応時間 (min)

3.05: 反応総液量 (ml)

0.05: 反応に供した酵素試料液量 (ml)

X: 酵素試料液の検品濃度 (mg/ml)